

# **Thetis** **Medical**


## **Achilles Splint Specialists**

Thetis splints improve the recovery experience for patients, save hospital staff time and reduce costs.



# Achilles Trauma Splint

## Intended Use

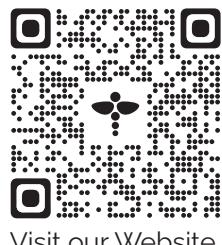

Visit our Website

**Hospitals** - The splint is applied in A&E/ER. Application requires no training. The Thetis Trauma Splint holds the foot plantar flexed until the patient is seen by a specialist.

**Physiotherapists** - Provides sports team medical professionals with easy pitch-side splintage allowing direct referral to specialist & imaging. This way the patient avoids a long and unnecessary A&E wait.

## Problem

Currently, patients who visit A&E with a suspected Achilles rupture are placed in a cast. This cast needs to be removed (and sometimes reapplied) if the patient requires a scan or operation.

- ✖ 5 Layer application
- ✖ Staff require training
- ✖ Difficult to apply
- ✖ 15 Minute drying time

## Solution

By using a comfortable, padded, plastic shell wrapped in crepe bandage, the process is transformed.

- ✓ 2 Layer application
- ✓ No training required
- ✓ Easier to apply
- ✓ Reduce patient waiting time



The Thetis Trauma Splint saves patients and staff time in A&E. It also facilitates an efficient pathway for scanning and making decisions about treatment


Musculo-Skeletal Radiologist, Past UKRCO Vice President and BSSR Council

# Application


SINGLE PATIENT USE ONLY.

BEFORE YOU START YOU WILL NEED:

1. THETIS SPLINT - Choose Right or Left, Small or Large
2. Crepe Bandage, Tape and Scissors (or Cohesive Bandage).

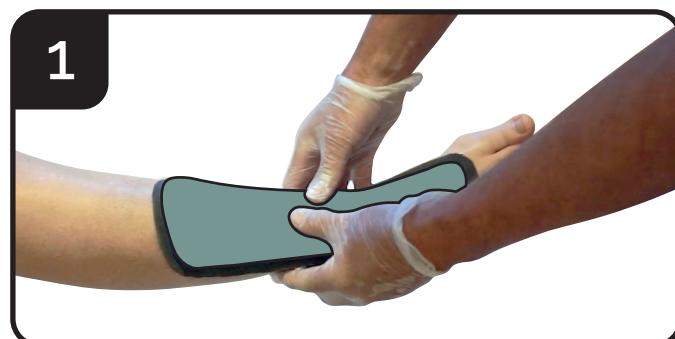

1. Ask the patient to point their toes down (tip-toes). Position the splint on the lower shin so that it also covers the top of the foot.

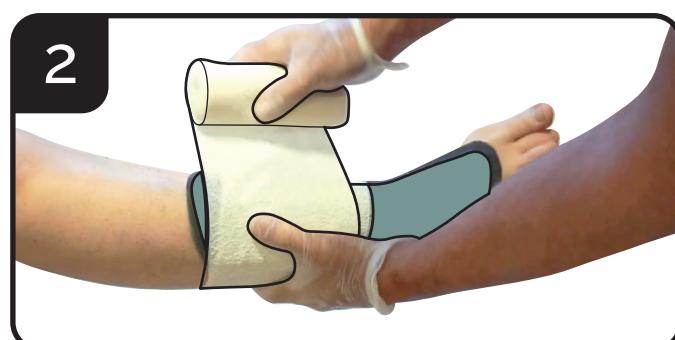


2. The foam of the lower end of the splint should sit just above the joints where the toes meet the foot. Wrap with crepe bandage firmly but not too tight.



3. Fix the bandage with medical tape.

## Splint Sizing

	<b>SMALL</b>	<b>LARGE</b>
UK	< 8	≥ 8
EU	< 42	≥ 42
US men	< 9	≥ 9
US women	< 10.5	≥ 10.5

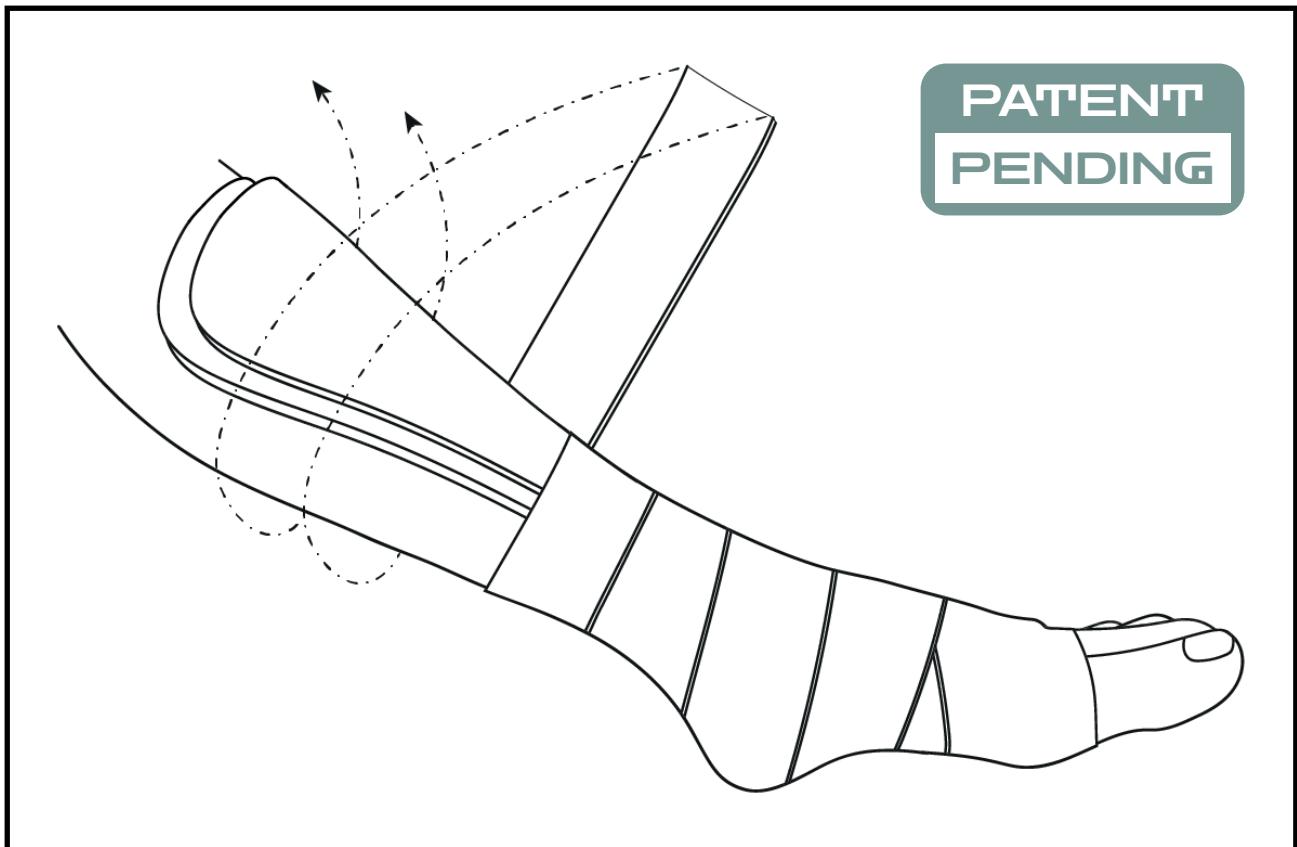


---

## Quality

---

Thetis Medical Ltd is a UK based company. All our Achilles splints are made domestically. Our manufacturers are ISO 9001 certified.

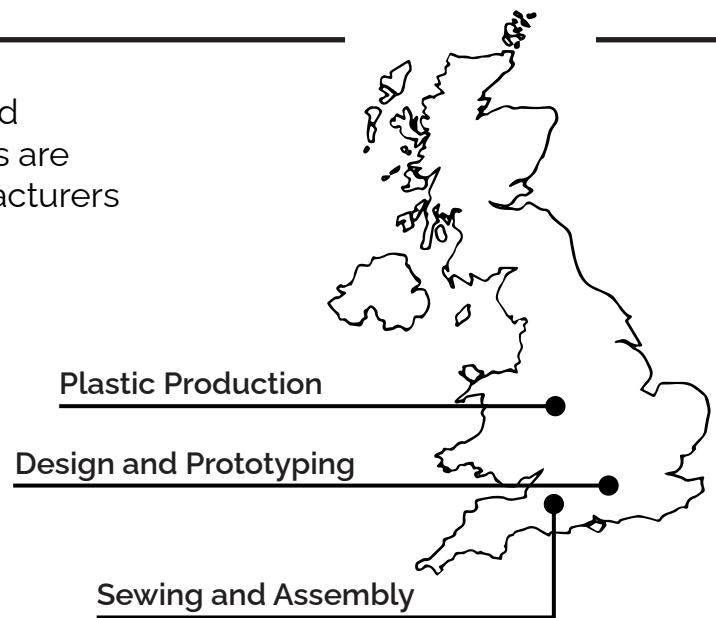

**UK  
CA**

[info@thetismedical.com](mailto:info@thetismedical.com)

