

Thetis **Medical**


Achilles Splint Specialists

Thetis splints improve the recovery experience for patients, save hospital staff time and reduce costs.



Achilles Trauma Splint

Intended Use

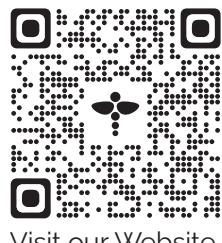

Visit our Website

Hospitals - The splint is applied in A&E/ER. Application requires no training. The Thetis Trauma Splint holds the foot plantar flexed until the patient is seen by a specialist.

Physiotherapists - Provides sports team medical professionals with easy pitch-side splintage allowing direct referral to specialist & imaging. This way the patient avoids a long and unnecessary A&E wait.

Problem

Currently, patients who visit A&E with a suspected Achilles rupture are placed in a cast. This cast needs to be removed (and sometimes reapplied) if the patient requires a scan or operation.

- ✖ 5 Layer application
- ✖ Staff require training
- ✖ Difficult to apply
- ✖ 15 Minute drying time

Solution

By using a comfortable, padded, plastic shell wrapped in crepe bandage, the process is transformed.

- ✓ 2 Layer application
- ✓ No training required
- ✓ Easier to apply
- ✓ Reduce patient waiting time



The Thetis Trauma Splint saves patients and staff time in A&E. It also facilitates an efficient pathway for scanning and making decisions about treatment

Musculo-Skeletal Radiologist, Past UKRCO Vice President and BSSR Council

Application


BEFORE YOU START YOU WILL NEED:

1. THETIS SPLINT - Choose Right or Left, Small or Large
2. Crepe Bandage, Tape and Scissors (or Cohesive Bandage).




Video Demonstration

1

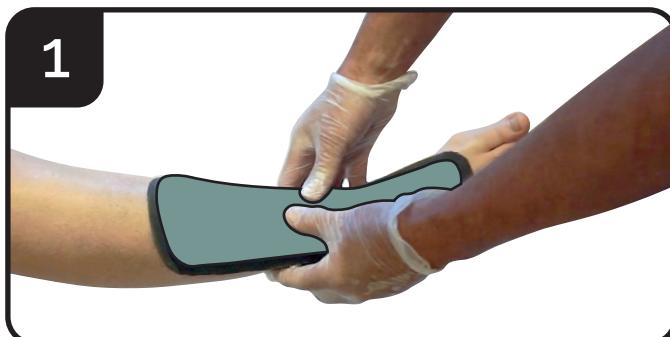

1. Ask the patient to point their toes down (tip-toes). Position the splint on the lower shin so that it also covers the top of the foot.

2

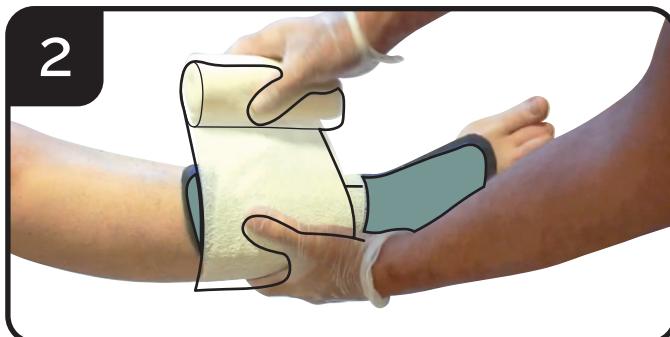

2. The foam of the lower end of the splint should sit just above the joints where the toes meet the foot. Wrap with crepe bandage firmly but not too tight.


3



3. Fix the bandage with medical tape.


Splint Sizing

	SMALL	LARGE
UK	< 8	≥ 8
EU	< 42	≥ 42
US men	< 9	≥ 9
US women	< 10.5	≥ 10.5

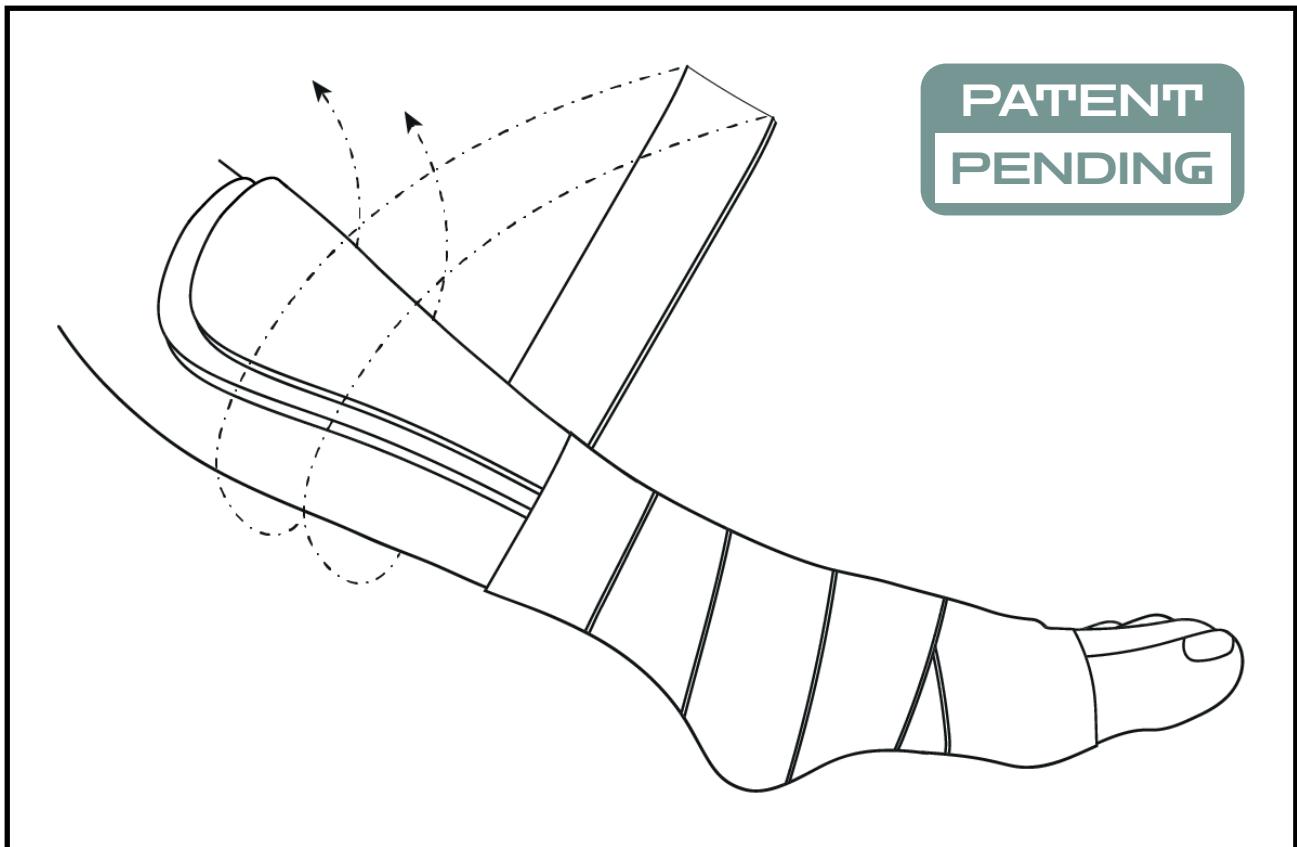

Quality

Thetis Medical Ltd is a UK based company. All our Achilles splints are made domestically. Our manufacturers are ISO 9001 certified.

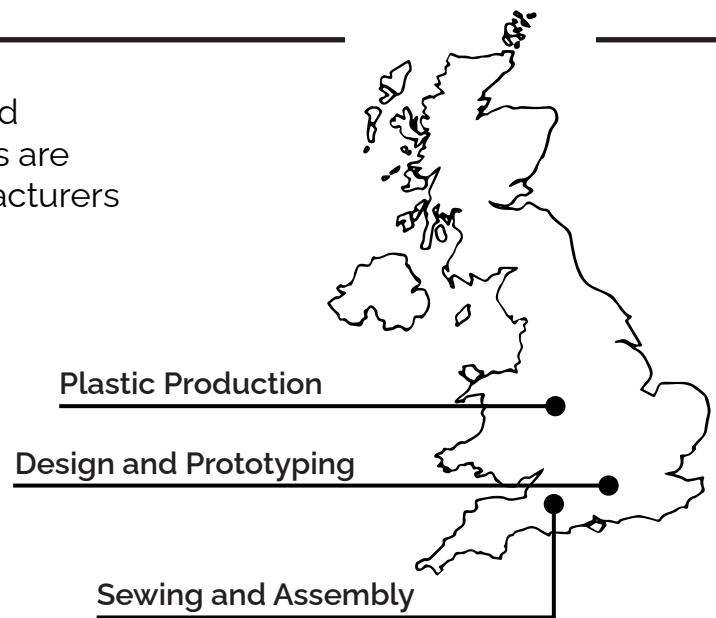

**UK
CA**

info@thetismedical.com

